

# Kansas Agricultural Experiment Station Research Reports

---

Volume 0

Issue 1 *Cattleman's Day* (1993-2014)

Article 76

---

2014

## Evaluating the effectiveness of transport media on Shiga toxin-producing *E. coli* serotypes

Nicholas W. Baumann

A. West

Randall K. Phebus

Follow this and additional works at: <https://newprairiepress.org/kaesrr>

 Part of the [Other Animal Sciences Commons](#)

---

### Recommended Citation

Baumann, Nicholas W.; West, A.; and Phebus, Randall K. (2014) "Evaluating the effectiveness of transport media on Shiga toxin-producing *E. coli* serotypes," *Kansas Agricultural Experiment Station Research Reports*: Vol. 0: Iss. 1. <https://doi.org/10.4148/2378-5977.1479>

This report is brought to you for free and open access by New Prairie Press. It has been accepted for inclusion in Kansas Agricultural Experiment Station Research Reports by an authorized administrator of New Prairie Press. Copyright 2014 Kansas State University Agricultural Experiment Station and Cooperative Extension Service. Contents of this publication may be freely reproduced for educational purposes. All other rights reserved. Brand names appearing in this publication are for product identification purposes only. No endorsement is intended, nor is criticism implied of similar products not mentioned. K-State Research and Extension is an equal opportunity provider and employer.



---

# Evaluating the effectiveness of transport media on Shiga toxin-producing *E. coli* serotypes

## Abstract

One of the key issues involved in accurately testing beef and the environment for the presence of specific bacteria, particularly pathogens such as Shiga toxin-producing *Escherichia coli* (STEC), is maintaining the viability of the microorganisms when transporting samples from the field to the laboratory. This process may take up to three days when considering collection, shipping and laboratory preparation times. Allowing the target bacteria to increase or decrease in numbers during transit is undesirable, so samples must be kept chilled and the media used for transport must offer a stable but non-nutritive environment. Three commonly used non-selective transport media were evaluated for their ability to maintain original STEC levels during transport. Holding temperature may vary during shipping, so this study evaluated two separate temperatures as co-variables.

## Keywords

Cattlemen's Day, 2014; Kansas Agricultural Experiment Station contribution; no. 14-262-S; Report of progress (Kansas State University. Agricultural Experiment Station and Cooperative Extension Service); 1101; Beef Cattle Research, 2014 is known as Cattlemen's Day, 2014; Beef; Shiga toxin-producing *E. coli*; Transportation; Temperature

## Creative Commons License



This work is licensed under a [Creative Commons Attribution 4.0 License](https://creativecommons.org/licenses/by/4.0/).

# Evaluating the Effectiveness of Transport Media on Shiga Toxin-Producing *E. coli* Serotypes

*N. Baumann, A. West, and R.K. Phebus*

## Introduction

One of the key issues involved in accurately testing beef and the environment for the presence of specific bacteria, particularly pathogens such as Shiga toxin-producing *Escherichia coli* (STEC), is maintaining the viability of the microorganisms when transporting samples from the field to the laboratory. This process may take up to three days when considering collection, shipping and laboratory preparation times. Allowing the target bacteria to increase or decrease in numbers during transit is undesirable, so samples must be kept chilled and the media used for transport must offer a stable but non-nutritive environment. Three commonly used non-selective transport media were evaluated for their ability to maintain original STEC levels during transport. Holding temperature may vary during shipping, so this study evaluated two separate temperatures as co-variables.

## Experimental Procedures

Buffered peptone water (BPW; Becton, Dickinson and Company, Sparks, MD), Cary-Blair transport media (CB; Oxoid LTD, Basingstoke, Hampshire, England), and maximum recovery diluent (MR; Oxoid LTD) media were selected for evaluation. Sets of these media types were individually inoculated with STEC serotypes O26, O45, O103, O111, O145, and O157 before being placed in refrigerators at 30 °F or 50 °F. Samples were removed and plated from each media/strain/temperature combination at 0, 12, 24, 48, and 72 hours. Plates were placed in an incubator at 99 °F for 24 hours before enumeration and comparison.

## Results and Discussion

The average log counts of each serotype on each media type at 39 °F and 50 °F are shown in Figures 1–6.

When transport medium samples were held at the desired refrigeration temperature 39 °F, all media types performed similarly; at the simulated abuse temperature of 50 °F, the Cary-Blair medium was the only one that maintained the bacteria at close to the original inoculation level.

## Implications

The use of Cary-Blair transport medium in conjunction with proper chilled storage temperatures will help ensure that accurate results are obtained for field-based research studies, and for daily beef processing samples being transported to analytical laboratories for enumeration of STEC.

MEAT AND FOOD SAFETY

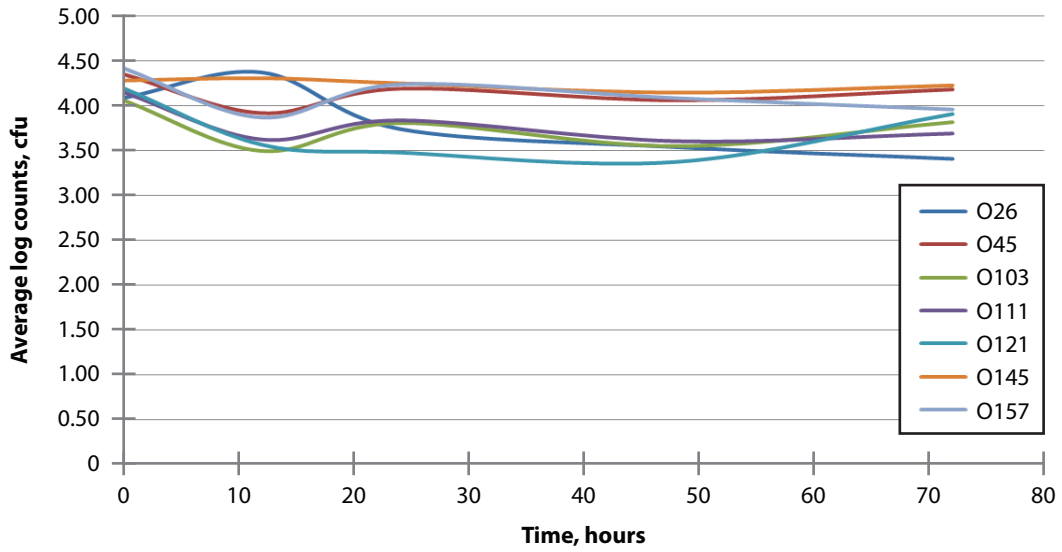


Figure 1. The average log count (cfu) of E. coli O26, O45, O103, O111, O145, and O157 on buffered peptone water and incubated at 39°F for up to 72 hours.

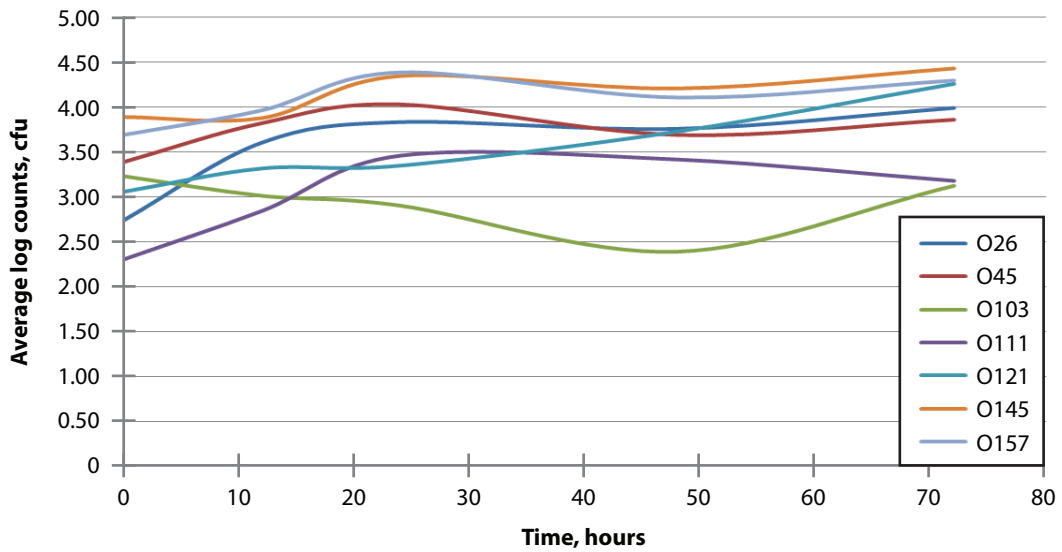


Figure 2. The average log count (cfu) of E. coli O26, O45, O103, O111, O145, and O157 on Cary-Blair transport media and incubated at 39°F for up to 72 hours.

MEAT AND FOOD SAFETY

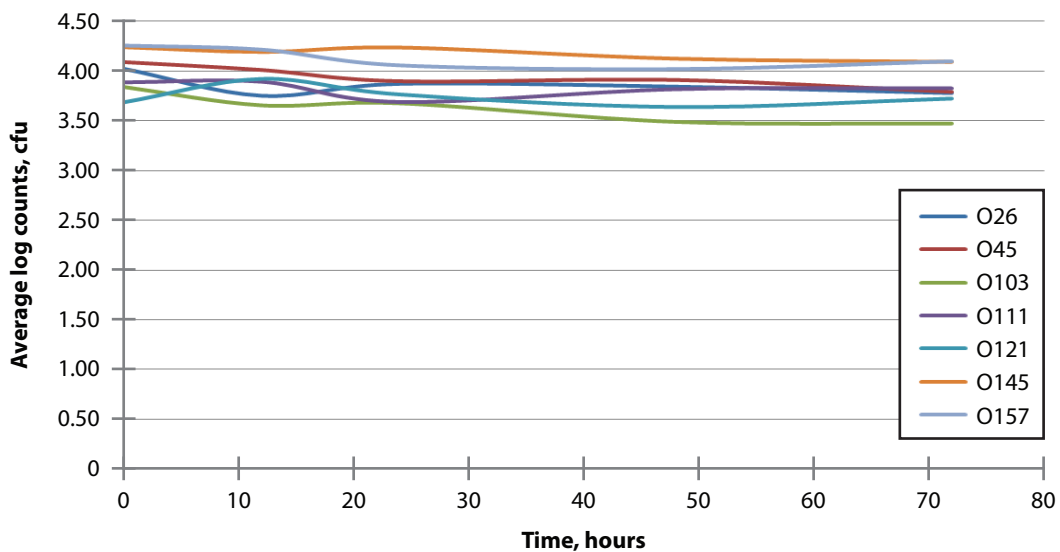


Figure 3. The average log count (cfu) of E. coli O26, O45, O103, O111, O145, and O157 on maximum recovery diluent and incubated at 39° F for up to 72 hours.

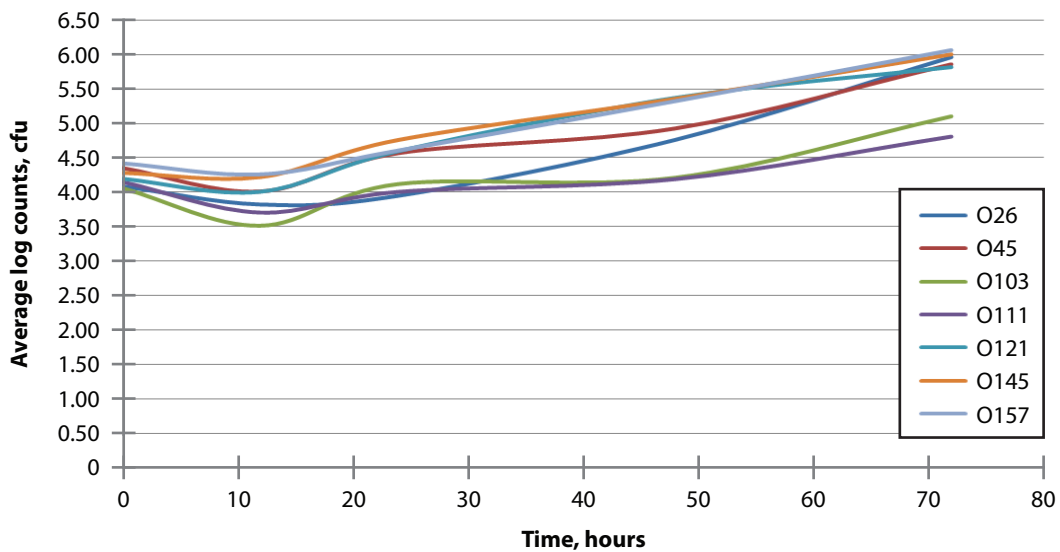


Figure 4. The average log count (cfu) of E. coli O26, O45, O103, O111, O145, and O157 on buffered peptone water and incubated at 50° F for up to 72 hours.

## MEAT AND FOOD SAFETY

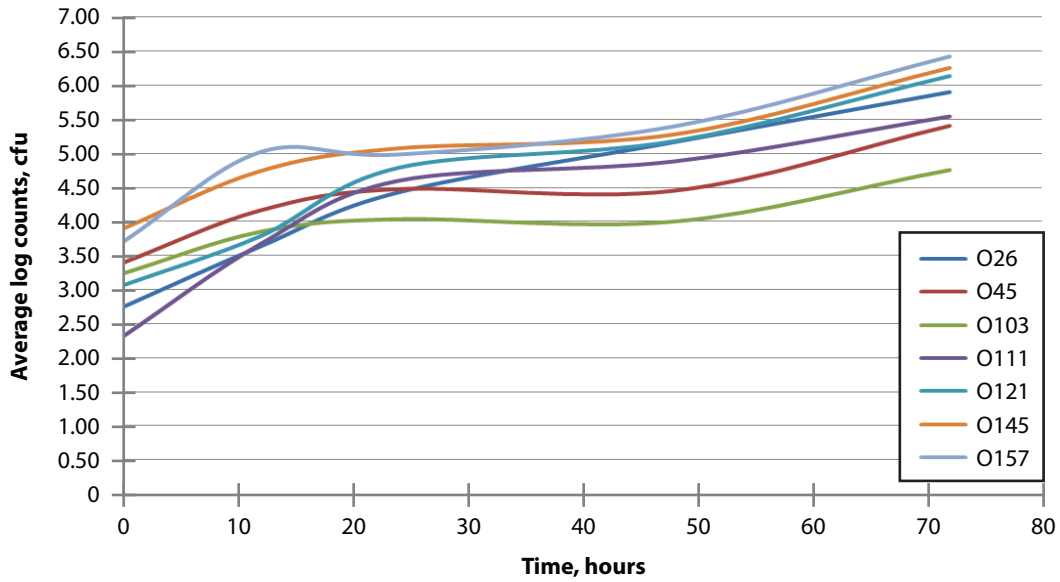


Figure 5. The average log count (cfu) of E. coli O26, O45, O103, O111, O145, and O157 on Cary-Blair transport media and incubated at 50° F for up to 72 hours.

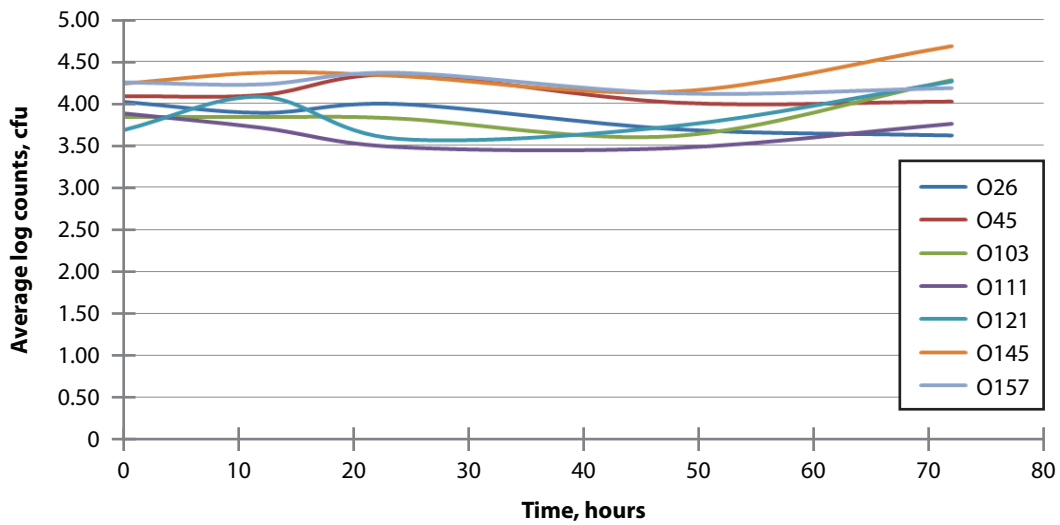


Figure 6. The average log count (cfu) of E. coli O26, O45, O103, O111, O145, and O157 on maximum recovery diluent and incubated at 50° F for up to 72 hours.